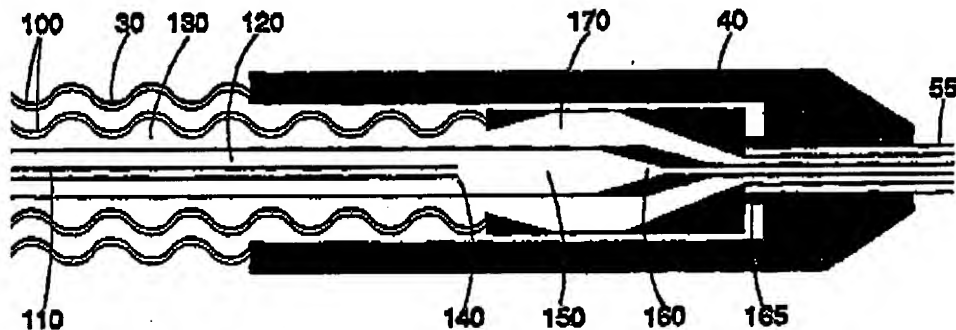


PCTWORLD INTELLECTUAL PROPERTY ORGANIZATION
International Bureau

INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification ⁶ : A61B 17/36		A2	(11) International Publication Number: WO 95/30379
			(43) International Publication Date: 16 November 1995 (16.11.95)
(21) International Application Number: PCT/GB95/00739		(81) Designated States: JP, US, European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE).	
(22) International Filing Date: 31 March 1995 (31.03.95)		Published <i>Without international search report and to be republished upon receipt of that report.</i>	
(30) Priority Data: 9409229.3 10 May 1994 (10.05.94) GB 9409230.1 10 May 1994 (10.05.94) GB 9409231.9 10 May 1994 (10.05.94) GB			
(71) Applicant (for all designated States except US): SPEMBLY MEDICAL LIMITED [GB/GB]; Newbury Road, Andover, Hampshire SP10 4DR (GB).			
(72) Inventors; and (75) Inventors/Applicants (for US only): GODDARD, Robert, William [GB/GB]; Kalmina, 74B The Crescent, Andover, Hampshire SP10 3BU (GB). PENNY, Leslie, John [GB/GB]; 215 Weyhill Road, Andover, Hampshire SP10 3LJ (GB).			
(74) Agents: TURNER, James, Arthur; D. Young & Co., 21 New Fetter Lane, London EC4A 1DA (GB).			

(54) Title: **CRYOSURGICAL INSTRUMENT**

(57) Abstract

A cryosurgical instrument comprises a cooling tip operable to be cooled by a flow of a cryogenic fluid within the tip; and a fluid transfer chamber having a fluid inlet for receiving cryogenic fluid, a first fluid outlet communicating with the cooling tip for supplying cryogenic fluid to the cooling tip, and a second fluid outlet communicating with an atmospheric air vent.

BEST AVAILABLE COPY

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Mauritania
AU	Australia	GE	Georgia	MW	Malawi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BG	Bulgaria	IE	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Belarus	KE	Kenya	RO	Romania
CA	Canada	KG	Kyrgyzstan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic of Korea	SD	Sudan
CG	Congo	KR	Republic of Korea	SE	Sweden
CH	Switzerland	KZ	Kazakhstan	SI	Slovenia
CI	Côte d'Ivoire	LI	Liechtenstein	SK	Slovakia
CM	Cameroon	LK	Sri Lanka	SN	Senegal
CN	China	LU	Luxembourg	TD	Chad
CS	Czechoslovakia	LV	Latvia	TG	Togo
CZ	Czech Republic	MC	Monaco	TJ	Tajikistan
DE	Germany	MD	Republic of Moldova	TT	Trinidad and Tobago
DK	Denmark	MG	Madagascar	UA	Ukraine
ES	Spain	ML	Mali	US	United States of America
FI	Finland	MN	Mongolia	UZ	Uzbekistan
FR	France			VN	Viet Nam
GA	Gabon				

WO 93/30379

PCT/GB93/00739

1

CRYOSURGICAL INSTRUMENT

This invention relates to cryosurgical instruments.

It is known to use cryosurgical instruments in surgical operations in which
5 living tissue is destroyed (or ablated) by cooling the tissue to well below the freezing
point of water (e.g. -20° to -30° Centigrade)

A common form of a cryosurgical instrument for use in these surgical
procedures comprises a hand-held cryosurgical probe, connected by various cryogenic
fluid supply and exhaust conduits to a supply of cryogenic fluid such as a dewar
10 vessel. A cooling tip of the probe is inserted into or touched against the tissue to be
ablated and is then cooled by passing the cryogenic fluid into the interior of the
cooling tip.

Often, the cryosurgical probe comprises a handle portion (for the surgeon or
operator to hold) connected to the cooling tip by a relatively slender uncooled shaft.
15 The purpose of the slender shaft is to allow the surgeon, when necessary, to insert the
cooling tip deep into the tissue to be ablated or to manoeuvre the tip into a deep
recess within the patient's body. Despite its potentially narrow cross section, the shaft
carries at least a cryogenic fluid supply conduit for supplying cryogenic fluid to the
cooling tip and an exhaust conduit for carrying cryogenic fluid from the cooling tip.

20 However, the very fact that the shaft has to be slender can cause problems with
the flow of cryogenic fluid within the shaft, particularly along the supply conduit in
a direction towards the cooling tip. This problem arises because the supply conduit
has to be made narrow in order to fit into the shaft, and so a high cryogenic fluid
velocity is required to force the fluid along the supply conduit. However, when the
25 instrument is first used after a period of non-use, all of the conduits (including the
supply conduit in the shaft) will tend to be warm. This can lead to localised boiling
or expansion of the cryogenic fluid and it can be difficult to establish a flow rate of
the cryogenic fluid towards the tip which is sufficient to cool the whole system down
and to provide a high enough fluid velocity that the cryogenic fluid will pass along
30 the slender shaft to reach the cooling tip.

WO 93/04647 discloses a cryosurgical instrument in which a metal tip of a
cryosurgical probe is cooled by supplying liquid nitrogen to the interior of the tip.

WO 95/30379

PCT/GB95/00739

2

This instrument provides small vent holes to allow gas bubbles formed in the cryogenic fluid supply conduit to vent to an exhaust conduit before the bubbles reach the cooling tip.

5 This invention provides a cryosurgical instrument comprising: a cooling tip operable to be cooled by a flow of a cryogenic fluid within the tip; and a fluid transfer chamber having a fluid inlet for receiving cryogenic fluid, a first fluid outlet communicating with the cooling tip for supplying cryogenic fluid to the cooling tip, and a second fluid outlet communicating with an atmospheric air vent.

10 The invention addresses the problem of establish a sufficient flow rate of the cryogenic fluid towards the tip, particularly when the instrument is first used, by providing two possible return paths for the cryogenic fluid. Cryogenic fluid reaching the fluid transfer chamber can either pass through the first fluid outlet towards the cooling tip or can return to the atmospheric air vent via the second fluid outlet (which is preferably separate from an exhaust passage for fluid returning from cooling
15 operation in the tip). This extra return path means that more fluid is passing along the various supply conduits towards the fluid transfer chamber, and so the fluid velocity is higher. This increased fluid velocity or flow rate can assist in forcing the fluid along a narrow conduit towards the cooling tip (if necessary) and can also help to cool the whole system down when the system is first used.

20 Preferably the instrument comprises a handle portion housing the fluid transfer chamber; and an elongate shaft linking the handle portion and the cooling tip. By housing the fluid transfer chamber in the handle portion, the extra fluid flow produced by this arrangement passes along all of the supply conduits leading to the handle portion.

25 In order to improve the proportion of the cryogenic fluid entering the fluid transfer chamber which then passes through the first fluid outlet, it is preferred that the fluid inlet and the first fluid outlet are disposed at opposite sides of the fluid transfer chamber. This proportion can further be improved in a preferred embodiment in which the fluid inlet is positioned to direct cryogenic fluid towards the first fluid
30 outlet. In order to direct the fluid into the first fluid outlet, it is preferred that the interior of the fluid transfer chamber is tapered towards the first fluid outlet.

Preferably the fluid inlet and the second fluid outlet are disposed on the same

WO 95/30379

PCT/GB95/00739

3

side of the fluid transfer chamber. This can encourage gas bubbles formed in the supply conduit to vent via the second outlet, because the liquid's greater momentum will tend to cause the liquid to cross the chamber towards the first outlet. This can increase the proportion of liquid at the cooling tip, thus improving the cooling efficiency.

In order that fluid returning via the second fluid outlet to the atmospheric air vent can be used efficiently to cool the supply conduit carrying the fluid to the fluid transfer chamber, it is preferred that the instrument comprises an elongate supply conduit for supplying cryogenic fluid to the fluid inlet; and a by-pass return conduit substantially coaxial with the supply conduit for carrying cryogenic fluid from the second fluid outlet to the atmospheric air vent.

In order to avoid an unnecessary waste of cryogenic fluid when an appropriate flow rate has been established, the instrument preferably comprises a valve for controlling the flow of fluid through the atmospheric air vent. (In contrast, the apparatus of WO 93/04647 would in practice require a liquid scavenging system to recover cryogenic liquid exhausted through the vent holes).

Preferably the instrument comprises means for opening the valve to allow the flow of fluid through the atmospheric air vent for an initial predetermined period of use of the instrument.

It is preferred that the predetermined period is between about three minutes and about five minutes. After this time the flow of fluid through the second outlet can be completely or substantially stopped by the valve.

Since problems involving an insufficient fluid flow rate might arise after the predetermined period has elapsed, it is preferred that the instrument comprises means for detecting a flow rate of cryogenic fluid through the instrument; and means for opening the valve to allow the flow of fluid through the atmospheric air vent if the flow rate is less than a predetermined amount.

In order to avoid the potential hazard of liquid cryogenic fluid being vented directly to the atmosphere, it is preferred that the atmospheric air vent comprises a vented liquid trap chamber.

If any of the conduits used to transport the supply and exhaust cryogenic fluid become contaminated with moisture or debris, the performance of the instrument can

WO 95/30379

PCT/GB95/00739

4

suffer. In the case of particulate debris, problems can arise if the flow of cryogenic fluid to the tip is impeded. However, moisture ingress is a more serious problem, since any moisture present in the conduits when the instrument is first used can freeze and potentially block one of the conduits.

5 Accordingly, in a preferred embodiment, the instrument has an exhaust conduit for carrying cryogenic fluid from the tip to an exhaust atmospheric air vent; and means for directing a flow of substantially dry purging gas along at least a part of the exhaust conduit to exhaust through the atmospheric air vent.

10 This feature addresses the potential problem of moisture ingress into the system when cooling of the tip is not being performed, in particular through the exhaust vent, by maintaining a small but steady flow of purging gas out of the exhaust vent. This flow can be maintained when the system is not in use. In fact, preferably, the flow is only maintained when the system is not in use.

15 In order to avoid moisture ingress into the potentially very narrow fluid passages in the cooling tip, it is preferred that the directing means is also operable to direct the flow of gas through the cooling tip. This can also be advantageous in a reusable probe, since some moisture can become trapped in the cooling tip during a steam sterilisation process. The purging gas can help to drive out this trapped moisture.

20 Preferably the directing means is operable to direct the flow of gas into the by-pass fluid return conduit (i.e. towards the second fluid outlet in the fluid transfer chamber, when the probe is not in cooling operation). This arrangement allows the fluid supply conduit to be sealed when the instrument is not in use, but still allowing the flow of purging gas to the interior of the tip.

25 The gas and the cryogenic fluid can conveniently be generated from the same source. This can give an added advantage that the steady flow of the purging gas can have a cooling effect on parts of the fluid supply such as a dip tube in a dewar vessel. It is therefore preferred that the instrument comprises a dewar vessel for storing a liquefied cryogenic fluid; and a dip tube extending into the fluid in the dewar vessel,
30 the dip tube having a first passage for carrying cryogenic fluid for use in cooling the cooling tip, and a second passage for carrying gaseous cryogenic fluid into the third conduit. Since the purging gas is generated by gently boiling off liquid cryogenic

WO 95/30379

PCT/GB95/00739

5

fluid when the dewar is pressurised, it is very dry indeed. It has been found that such extremely dry gas can actually tend to absorb moisture in the tubes.

In order to regulate the flow of the purging gas through the system, it is preferred that the instrument comprises a flow restricting valve in the path of the gas flow.

As an added safety feature, it is preferred that the instrument comprises means for detecting the pressure of the purging gas. Preferably the instrument also comprises means for inhibiting cooling operation of the instrument if the detected pressure is outside a predetermined range of gas pressures.

These safety features rely on detecting changes in pressure in the by-pass passages. When the system is not in use, with purging gas flowing along the by-pass passages, if the probe is not connected the gas pressure along those passages will be substantially atmospheric (10^5 Pa). When a correctly functioning probe is connected into the system, the gas pressure detected by the pressure detector will rise. If the probe is connected but has a blockage, the detected pressure will rise further.

Accordingly, the output of the pressure detector can be compared with a first threshold (e.g. 10^4 Pa over atmospheric pressure, or alternatively 1.1×10^5 Pa), with cooling operation being inhibited if the detected pressure is below that threshold. A warning light could be illuminated on the control console to indicate that the probe is not connected. (This could be particularly useful in a multi-channel system where a probe could be connected to one of several fluid outlets; the system would prevent fluid being sent to an outlet when the probe was accidentally connected to another outlet).

Similarly, if the detected pressure exceeds another threshold (e.g. 2×10^4 Pa over atmospheric pressure or alternatively 1.2×10^4 Pa) then operation could again be inhibited, with a further warning light indicating a probe fault being illuminated.

Preferably the cryogenic fluid is liquid nitrogen and the gas is nitrogen gas.

When the instrument is used in a surgical procedure, the instrument has to be kept sterile to avoid the risk of infecting the patient during the procedure.

In use, the instrument has to be connected to a supply source of cryogenic fluid. It would not be feasible to sterilise the entire supply source and instrument combination, so the instrument is generally fitted with a connector to allow temporary

WO 95/30379

PCT/GB95/00739

6

connection of the instrument to the supply source during a surgical procedure. The instrument (and possibly the connector) can then be sterilised before use.

5 A preferred embodiment of the invention provides connection apparatus comprising: complementary first and second connectors having cooperable locking means operable to lock the first connector to the second connector; in which: the first connector comprises means for activating the locking means to lock the first connector to the second connector; and the second connector comprises means for releasing the locking means to release the first connector from the second connector.

10 This feature recognises that two potentially conflicting requirements apply to connectors for use in sterile environments. First, the connector should be lockable to avoid the potentially hazardous possibility of the connector becoming detached from the supply source during the surgical procedure. However, if the locking mechanism requires a sterile operator or scrub nurse to touch the 'supply' (non-sterile) side of the connection, then that operator's hands can become contaminated with non-sterile
15 contaminants. This means that the operator then has to go through a lengthy washing process and, in many cases, put on new sterile clothing and gloves before the operator can take a further part in the surgical procedure.

20 These apparently conflicting requirements are addressed by providing a connection apparatus in which two connectors can be locked together by a sterile operator holding only one of the connectors (the sterile connector). For security, the release mechanism is not mounted on the sterile connector. When it is time to release the connectors, the surgical procedure has been completed so a non-sterile operator can use the release means on the other (non-sterile) connector to release the locking of the connectors without concern over contamination with non-sterile material.

25 If the probe needs to be disconnected during a sterile procedure, for example to move the probe from one fluid supply to another in a multi-channel system, a sterile operator needs to operate the sterile side of the connector and a non-sterile operator needs to release the other side of the connector. This dual-operator interlock can help prevent inappropriate disconnection of the probe during a procedure.

30 Preferably the first and second connectors comprise complementary bayonet locking formations.

In one preferred embodiment the second connector comprises a sleeve

WO 95/30379

PCT/GB95/00739

7

resiliently biased to overlie and prevent release of the bayonet locking formations, the first connector having sleeve displacing formations to displace the sleeve when the first and second connectors are locked together. In this way, the action of pushing the two connectors together displaces the sleeve to allow the bayonet to be locked. The sleeve then overlies the bayonet formations and prevents the two connectors being unlocked until the sleeve is moved by an operator.

In order to reduce further the risk of operator contamination, it is preferred that the first connector comprises guard formations for at least partially shielding an operator's hand when the first and second connectors are locked together.

Because the connectors are used to carry cryogenic fluids, it is preferred that the apparatus comprises a heater for heating at least one of the first connector and the second connector. The heater reduces the chance of moisture in or on the connector freezing and temporarily preventing the connector from being disconnected.

Viewed from a second aspect this invention provides a cryosurgical instrument comprising: a probe having a cooling tip operable to be cooled by a flow of a cryogenic fluid within the tip; means for providing a return path for cryogenic fluid from the probe without the fluid passing into the cooling tip; and means for preventing fluid flow via the return path after a predetermined period of operation of the instrument.

The invention will now be described by way of example with reference to the accompanying drawings, throughout which like parts are referred to by like references, and in which:

Figure 1 is a schematic diagram of a cryosurgical probe system;

Figure 2 is a schematic diagram of a cryosurgical probe handle;

Figures 3a and 3b are enlarged schematic diagrams of a cryosurgical probe tip;

Figures 4a, 4b and 4c are schematic diagrams of a cryogenic fluid connector;

Figures 5a and 5b are schematic diagrams illustrating the operation of the cryogenic fluid connector of Figures 4a, 4b and 4c; and

Figure 6 is a schematic diagram of a cryogenic fluid supply apparatus.

Figure 1 is a schematic diagram of a cryosurgical probe system. The cryosurgical probe system is suitable for use in surgical operations in which living tissue is destroyed (or ablated) by cooling the tissue to well below the freezing point

WO 95/30379

PCT/GB95/00739

8

of water (e.g. -20° to -30° Centigrade).

A cryogenic fluid, such as liquid nitrogen, is stored in a dewar vessel 10 forming part of a fluid supply apparatus 5. Alternative cryogenic fluids include freon 14, freon 22, freon 13, liquefied air and normally gaseous hydrocarbons (although inert fluids are preferred to flammable fluids).

The cryogenic fluid is supplied via an overhead connecting tube 20 and a flexible connecting tube 30 to a probe handle 40 and eventually to a probe tip 50. The flexible connecting tube 30 is linked to the overhead connecting tube 20 by a releasable connector (see Figures 4a, 4b, 4c, 5a and 5b). This allows the probe and flexible connecting tube to be steam sterilised for re-use.

The cryogenic supply apparatus 5 also has a by-pass exhaust chamber 15. This is a liquid trapping chamber vented to the atmosphere. Its role is to prevent the potentially hazardous direct venting of liquid nitrogen from a by-pass return path (described below) to the atmosphere. Entry into the chamber 15 is controlled by a by-pass valve 25 mounted on the chamber 15.

When the probe system is in use during a surgical operation in an operating theatre, the fluid supply apparatus 5 is placed a short distance away from an operating table or other support for the patient. The overhead connecting tube 20 allows the transfer of cryogenic fluid between the fluid supply apparatus 5 and the cryosurgical probe held by the surgeon, without obstructing the surgeon's access to the patient. The overhead connecting tube 20 may be a rigid structure in its own right, or may be a flexible or semi-flexible tube supported on a suitable supporting frame (not shown).

In use, the probe tip is inserted into bodily tissue to be destroyed (for example, a part of the patient's liver) and cryogenic fluid is supplied from the dewar vessel 10 to the interior of the probe tip 50. This causes the probe tip 50 to be cooled, which in turn cools the bodily tissue into which the probe tip 50 is inserted. Over a period of a few minutes, a volume of the bodily tissue surrounding the probe tip 50 is cooled to about -20° to -30° Centigrade, killing the cells in that volume of tissue. When the required volume of tissue has been cooled in this way, the probe tip 50 can be allowed to thaw, either through conduction of bodily heat from the patient or by passing a heated gas into the probe tip 50. The heated gas can be generated by passing gas from a gas storage cylinder (not shown) through an electric heating element (not

WO 95/30379

PCT/GB95/00739

9

shown). The probe tip is then withdrawn from the bodily tissue. If necessary, the procedure outlined above may be repeated a number of times at different locations in the body or in a particular organ if a particularly large volume of tissue is to be destroyed.

5 The cryogenic fluid remains inside the probe tip 50 during the cooling process. In other words there is no potentially dangerous contact between the cryogenic fluid and the patient's bodily tissue.

Figure 2 is a schematic diagram of the probe handle 40, and Figures 3a and 3b are enlarged schematic diagrams of the cryosurgical probe tip 50.

10 Referring to Figures 2 and 3a, the probe handle 40 is connected to the flexible connecting tube 30 for the supply and exhaust of cryogenic fluid. The flexible connecting tube 30 is in fact a coaxial structure of several tubes. At the outside of the flexible connecting tube 30 are two skins 100, each formed of a double-walled structure of flexible corrugated tubing. This double-skinned structure provides good
15 heat insulation to prevent the outside of the flexible connecting tube 30 becoming hazarously cold when the probe is in use. Inside the inner skin 100 are two coaxial PTFE (polytetrafluoroethane) tubes, namely a cryogenic fluid supply tube 110 and, surrounding the cryogenic supply tube 110, a by-pass tube 120. A region 130 between the outside of the by-pass tube 120 and the inner skin 100 forms an exhaust
20 gas passage. Seals 165 are provided to terminate the annular space between the skins 100.

Cryogenic fluid supplied through the fluid supply tube 110 emerges at high velocity through a fluid outlet 140. The fluid emerges into a transfer chamber 150, from which some of the fluid may return towards the fluid supply apparatus along the
25 by-pass tube 120, and some may pass through a convergent nozzle 160 in a direction towards the probe tip 50.

For the first few minutes of operation of the probe system, the by-pass tube 120 is connected via the by-pass valve 25 (to be described below) to an atmospheric air vent. This means that a significant proportion of the cryogenic liquid emerging
30 from the fluid outlet 140 is able to return along the by-pass tube 120.

This arrangement means that the flow rate or fluid velocity of cryogenic fluid emerging from the outlet 140 can be increased by allowing some of the fluid to vent

WO 95/30379

PCT/GB95/00739

10

via the by-pass tube 120. The need for this arises because the fluid passing through the convergent nozzle 160 has to pass along a narrow fluid supply passage towards the tip 50. If the fluid velocity is increased by using the by-pass tube and valve, the initial flow of cryogenic fluid along the narrow fluid supply passage can be made much more reliable.

In addition, the by-pass system allows bubbles of gas formed in the fluid supply tube 110 to escape through the by-pass path, rather than passing the gas through the cooling tip 50. This diversion of the gas bubbles away from the tip can assist in cooling the tip down during the early stages of operation, since the cryogenic fluid has a much higher heat capacity per unit volume in liquid form than in gaseous form.

The separation of gas and liquid cryogenic fluid in the chamber 150 is helped by the relative positioning of the fluid outlet 140, the convergent nozzle 160 and the entrance to the by-pass tube 120. In a stream of gas and liquid cryogenic fluid, the liquid will have a much greater momentum than the gas, because it is about 600 times more dense than the gas phase. This means that the liquid is more likely to continue across the chamber 150 from the outlet 140 into the convergent nozzle 160, whereas the gas is more likely not to cross the chamber 150. The by-pass tube is connected to an atmospheric air vent and so provides a low-pressure, wide bore path for the gas to escape, whereas the convergent nozzle leads to a narrower bore, higher pressure path. Accordingly, the gas tends to follow the by-pass path.

The by-pass passages throughout the system are coaxial, surrounding the fluid supply tubes. This means that during the initial period of operation when the by-pass valve 25 is open, the fluid supply tubes are cooled by the flow of cryogenic fluid along the by-pass passages. Cooling the supply tubes is important to allow the cryogenic fluid to reach the probe tip 50 while still in liquid form (for efficient cooling of the probe tip 50), rather than as a gas.

After a predetermined time period in operation (e.g. 3 to 5 minutes) a control circuit (not shown) closes the by-pass valve 25 and so prevents further flow of cryogenic fluid along the by-pass tube 120. After this period in operation, the narrow fluid supply passage is full of the cryogenic fluid and is sufficiently cool that further fluid can flow easily down the passage. This means that the by-pass path is no longer

WO 95/30379

PCT/GB95/00739

11

needed, so it is shut off to avoid an unnecessary waste of the cryogenic fluid.

5 The by-pass valve 25 can also be controlled in response to the flow rate of exhaust fluid from the probe tip 50. If a problem arises which cuts down the flow of cryogenic fluid into the probe tip 50, this can be detected as a drop in the exhaust
10 flow rate below a threshold rate, by means of a suitable flow rate detector (not shown). In order to increase the fluid flow to the probe tip 50, the by-pass valve 25 can then be opened, either for a predetermined time or until the exhaust flow rate reaches the threshold rate. As described above, this has the effect of increasing the fluid velocity as the fluid approaches the convergent nozzle 160, and so increasing the
15 fluid flow rate into the probe tip 50.

During operation when the by-pass valve is closed, the consumption of cryosurgical fluid by the probe tip 50 is about 0.25 litres per minute.

Exhaust fluid emerging from the tip 50 passes into a coaxial exhaust chamber 170 which communicates with the exhaust passage 130.

15 Figure 3a is an enlarged schematic diagram of the probe tip 50. The probe tip 50 is connected to the handle 40 by a probe shaft 55. The probe shaft 55 has an outer wall 200 made of silver which surrounds a coaxial vacuum chamber 210. The vacuum chamber 210 provides a degree of heat insulation between the cryogenic fluid passing inside the shaft 55 and the bodily tissue in contact with the shaft 55. This
20 means that the freezing effect takes place primarily at the tip of the probe tip 50, rather than along the shaft 55, thus giving the surgeon better control over the location to be frozen. Within the shaft 55, there is a coaxial structure of a PTFE supply tube 220 and a surrounding exhaust passage 230. The returning exhaust fluid along the exhaust passage 230 can help to keep the incoming fluid cool.

25 The supply tube 220 continues as far as a heatsink 280 which then carries the cryogenic fluid to a chamber 240. The chamber 240 communicates with exhaust passages 250. Exhaust fluid passing through the exhaust passages 250 flows into the exhaust passage 230 through a series of spaced exhaust apertures 260. The path of fluid into the tip 50 and back along the exhaust passage 230 is shown schematically
30 in Figures 2 and 3a using direction arrows.

When the cryogenic fluid emerges into chamber 240, it starts to boil. There are two possible modes of boiling which can take place. If so-called 'film' boiling

WO 95/30379

PCT/GB95/00739

12

occurs, a film of nitrogen gas forms between the inner surface of the chamber 240 and the boiling liquid nitrogen. This film of gas reduces the heat conduction between the boiling liquid and the silver body of the tip 50 and thus reduces the cooling efficiency of the probe. The other possible mode of boiling is referred to as 'nucleate' boiling. In this mode boiling takes place at discrete bubble-forming sites at the liquid-solid boundary. The film of gas does not form, so the heat conduction between the boiling liquid and the outer silver body of the tip 50 is much better, and the cooling efficiency of the tip is thus improved.

It has been found that, unexpectedly, the nucleate boiling process can be encouraged by providing a thermally insulating coating on at least a part of the outer surface of the tip 50.

In the embodiment shown in Figure 3a, a thin PTFE coating 270 on the entire outer surface of the tip 50 is used. The coating is about 10 to 20×10^{-6} metres thick (preferably less than about 30×10^{-6} metres thick).

In an alternative embodiment shown schematically in Figure 3b, a similar insulating coating 275 is applied to a "boiling" region of the tip. This region is a part of the tip in which the liquid tends to boil, namely the passages along the milled flats of the screw-threaded connection 285. In this boiling region, the cryogenic liquid is in close thermal contact with tissue surrounding the tip, which contrasts with the insulated supply and exhaust passages 220, 230, where the liquid is relatively thermally isolated from ambient conditions. It has been found that coating the circumference of the tip over substantially this boiling region encourages the nucleate boiling referred to above.

Another factor which has been found to increase the likelihood of nucleate rather than film boiling is the construction of the cooling tip 50. Before the cryogenic fluid reaches the chamber 240, it passes inside the heatsink 280. (Before reaching the heatsink 280, the fluid was passing along the heat-insulating PTFE supply tube 220). The heatsink 280 is made of silver and is thermally connected by a screw thread connection (shown schematically as 285) to the outer silver casing 290 of the cooling tip. Longitudinal flats are milled on the heat sink 280 to form the exhaust passages 250. This technique of passing the cryogenic fluid either side of a heatsink thermally connected to the outer wall of the cooling tip 50 has been found to promote nucleate

WO 95/30379

PCT/GB95/00739

13

rather than film boiling within the cooling tip.

The outer silver wall 200 of the probe shaft 55 has a number of beneficial effects. It is in thermal contact with both the patient's tissue and the tip 50, and so it can conduct body heat towards the tip 50 to assist in defrosting the tip 50. It can act as a heat sink to help the shaft remain unfrozen even when the tip 50 is cooled. Also, the silver material helps to prevent corrosion of the shaft 55. This is important in a re-usable device.

Figures 4a, 4b and 4c are schematic diagrams of a cryogenic fluid connector. The connector is formed of two complementary connectors, a first connector 300 linked to the flexible connecting tube 30 and a second connector 310 linked to the overhead connecting tube 20. Figure 4a shows the first and second connectors mated together, Figure 4b shows the first connector and Figure 4c shows the second connector.

In Figure 4a, the path of fluid flow in the supply, by-pass and exhaust passages is indicated by arrows.

The connectors are locked together by a bayonet fixing formed by pins (350) on the first connector engaging in angled bayonet slots (not shown on Figures 4a, 4b and 4c) in the second connector. The two connectors are thus linked together by pushing the first connector into the second connector, and then twisting the first connector with respect to the second connector.

Sealing rings are provided at the mating surfaces 315 of the connectors.

The connector 310 has a resiliently biased sleeve 320 which overlies the bayonet connection and prevents the connection being disconnected. The operation of the sleeve 320 is illustrated in Figures 5a and 5b.

Figures 5a and 5b are schematic diagrams showing a bayonet slot 330 in the second connector, a retaining slot 340 formed by or connected to the sleeve 320, and a bayonet pin 350 of the first connector. (For clarity of Figures 4a to 4c, the retaining slot 340 has not been shown in these Figures. The retaining slot could be formed in the actual sleeve itself, or in a plate or other member attached to the sleeve).

When the pin 350 is pushed into the opening of the slot 330, it pushes against the spring-biased sleeve 320 until the pin has reached the bottom of the slot 330. The pin can then slide along the slot 330 during the rotational movement of the first and

WO 95/30379

PCT/GB95/00739

14

second connectors. During this movement the sleeve 320 is held back against the spring bias. However, when the pin reaches the retaining slot 340, the sleeve can return to its rest position and the pin is retained by the combination of the slot 330 and the retaining slot 340. The pin is no longer free to rotate, and so is captive
5 between the slot 330 and the slot 340. This situation is shown in Figure 5a.

Accordingly, to lock the connector as shown in Figure 5a, it is necessary for an operator to touch only the first connector 300. If this connector has been sterilised then a sterile operator or 'scrub nurse' can safely lock the two connectors together without having to touch the second (non-sterile) connector 310.

10 A flange 305 is provided to shield the operator's hand from accidentally touching the second connector during the connection process.

In order to unlock the two connectors, the sleeve 320 must be manually pushed against the spring bias into a position shown in Figure 5b. At this stage the surgical procedure is over, so it is safe for a (non-sterile) operator to touch the sleeve 320.
15 Once the sleeve has been pushed back, the first connector can be rotated and withdrawn to free the pin 350 from the slot 330. The sleeve is made of a plastics material, to reduce the danger of the sleeve freezing to the operator's hand when touched.

The connector 310 incorporates an electric heater (not shown), which is
20 energised at least during the time when the cryogenic fluid is passed through the connector. The heater reduces the chance of moisture in or on the connector freezing and temporarily preventing the connector from being disconnected.

Figure 6 is a schematic diagram of part of the cryogenic fluid supply apparatus
5 of Figure 1.

25 A dip tube 400 extends into the dewar vessel 10, to a level below the level of the cryogenic fluid in the dewar vessel. In normal operation, the cryogenic fluid to be supplied to the probe tip 50 passes along an axial supply passage 410, into the interior of a valve 420, and from the valve 420 into the axial supply passage 430 of the overhead connecting tube 20. The valve 420 can be opened and closed by means
30 of a solenoid controlled ruby or PTFE seal 425.

A second outlet from the dewar vessel is also provided, along an outer coaxial passage 440 in the dip tube 400. During quiescent operation when the dewar vessel

WO 95/30379

PCT/GB95/00739

15

is pressurised but the cooling tip is not being cooled (i.e. cryogenic liquid is not being supplied to the tip), this outlet provides a purging flow of cool gas through the various connecting tubes, the probe tip 50, and the exhaust passages of the cryosurgical probe system. The purging gas is actually generated by the liquid nitrogen in the dewar vessel boiling very gently.

The passage 440 is connected, via a flow restricting valve 445, to a chamber 450 which communicates with the by-pass return path 460 from the probe handle 40.

The by-pass valve 25 selectively vents or blocks fluid output from the chamber 450, as described above. When the probe system is not in use, the by-pass valve is closed and so purging gas flows at a low flow rate into the chamber 450 and then along the by-pass passages towards the probe handle 40 (i.e. in the opposite direction to the direction in which cryosurgical fluid flows along the by-pass passages when the by-pass valve is open). When the purging gas reaches the probe handle 40, it passes along the supply tube 220 of the probe shaft 55 to the probe tip 50 and then returns along the exhaust passages towards the fluid supply apparatus. The exhaust passages communicate with an exhaust vent 480 at the fluid supply apparatus.

The purging gas serves three main purposes. First, it can help to prevent the ingress of moisture into the system (in particular, through the exhaust vent 480), by maintaining a small but steady flow of gas out of the exhaust vent 480. Second, the small flow along the dip tube 400 can help to keep the dip tube cool, thus reducing the cool down time when the system is first used. Third, in a re-usable probe some moisture can become trapped in the cooling tip during a steam sterilisation process. The purging gas can help to drive out this trapped moisture. Since the purging gas is generated by boiling off liquid nitrogen, it is very dry indeed. It has been found that such extremely dry gas can actually tend to absorb moisture in the tubes.

A further use of the purging gas is obtained by placing a pressure sensor (not shown) to measure gas pressure along the by-pass passages, at any point between the gas outlet of the flow restricting valve 445 and the connector 310 (linked to the overhead connecting tube). For example, the pressure sensor could be connected to communicate with the chamber 450.

This pressure detector has two roles. It can be used to detect whether the probe is connected at the connector 310, as a safety feature to prevent cooling

WO 95/30379

PCT/GB95/00739

16

operation of the system being started without a probe in place to receive the cryogenic fluid. As an additional safety feature, the pressure detector can be used to detect blockages or malfunctions in the fluid passageways, before cooling operation is initiated.

5 These safety features rely on detecting changes in pressure in the by-pass passages. When the system is not in use, with purging gas flowing along the by-pass passages, if the probe is not connected the gas pressure along those passages will be substantially atmospheric (10^5 Pa). When a correctly functioning probe is connected into the system, the gas pressure detected by the pressure detector will rise. For
10 example, in a prototype of the present embodiment, the pressure in the by-pass passageways rises by about 100 millibars (10^4 Pa) when a correctly functioning probe is connected. If the probe is connected but has a blockage, the detected pressure will rise further.

 Accordingly, the output of the pressure detector can be compared with a first
15 threshold (e.g. 10^4 Pa over atmospheric pressure, or alternatively 1.1×10^5 Pa), with cooling operation being inhibited if the detected pressure is below that threshold. A warning light (not shown) could be illuminated on the control console to indicate that the probe is not connected. (This could be particularly useful in a multi-channel system where a probe could be connected to one of several fluid outlets; the system
20 would prevent fluid being sent to an outlet when the probe was accidentally connected to another outlet).

 Similarly, if the detected pressure exceeds another threshold (e.g. 2×10^4 Pa over atmospheric pressure or alternatively 1.2×10^4 Pa) then operation could again be inhibited, with a further warning light (not shown) indicating a probe fault being
25 illuminated.

WO 95/30379

PCT/GB95/00739
RECEIVED
CENTRAL FAX CENTER

17

NOV 08 2005

CLAIMS

1. A cryosurgical instrument comprising:
a cooling tip operable to be cooled by a flow of a cryogenic fluid within the
5 tip; and
a fluid transfer chamber having a fluid inlet for receiving cryogenic fluid, a
first fluid outlet communicating with the cooling tip for supplying cryogenic fluid to
the cooling tip, and a second fluid outlet communicating with an atmospheric air vent.
- 10 2. An instrument according to claim 1, comprising:
a handle portion housing the fluid transfer chamber; and
an elongate shaft linking the handle portion and the cooling tip.
- 15 3. An instrument according to claim 1 or claim 2, in which the fluid inlet and the
first fluid outlet are disposed at opposite sides of the fluid transfer chamber.
4. An instrument according to claim 3, in which the fluid inlet is positioned to
direct cryogenic fluid towards the first fluid outlet.
- 20 5. An instrument according to claim 4, in which the interior of the fluid transfer
chamber is tapered towards the first fluid outlet.
6. An instrument according to any one of the preceding claims, in which the fluid
inlet and the second fluid outlet are disposed on the same side of the fluid transfer
25 chamber.
7. An instrument according to any one of the preceding claims, comprising an
elongate supply conduit for supplying cryogenic fluid to the fluid inlet; and a return
conduit substantially coaxial with the supply conduit for carrying cryogenic fluid from
30 the second fluid outlet to the atmospheric air vent.
8. An instrument according to any one of the preceding claims, comprising a

WO 95/30379

PCT/GB95/00739

18

valve for controlling the flow of fluid through the atmospheric air vent.

9. An instrument according to claim 8, comprising means for opening the valve to allow the flow of fluid through the atmospheric air vent for an initial predetermined period of use of the instrument.

10. An instrument according to claim 9, in which the predetermined period is between about three minutes and about five minutes.

11. An instrument according to any one of claims 8 to 10, comprising means for detecting a flow rate of cryogenic fluid through the instrument; and means for opening the valve to allow the flow of fluid through the atmospheric air vent if the flow rate is less than a predetermined amount.

12. An instrument according to any one of the preceding claims, in which the atmospheric air vent comprises a vented liquid trap chamber.

13. A cryosurgical instrument according to any one of the preceding claims, comprising:

an exhaust conduit for carrying cryogenic fluid from the tip to an exhaust atmospheric air vent; and

means for directing a flow of substantially dry purging gas along at least a part of the exhaust conduit to exhaust through the exhaust atmospheric air vent.

14. An instrument according to claim 13, in which the directing means is operable to direct the flow of gas through the cooling tip.

15. An instrument according to claim 14, in which the directing means is operable to direct the flow of gas into the fluid return conduit.

16. An instrument according to any one of claims 13 to 15, comprising:
a dewar vessel for storing a liquefied cryogenic fluid; and

WO 95/30379

PCT/GB95/00739

19

a dip tube extending into the fluid in the dewar vessel, the dip tube having a first passage for carrying cryogenic fluid for use in cooling the cooling tip, and a second passage for carrying gaseous cryogenic fluid into the fluid return conduit.

5 17. An instrument according to claim 16, comprising a flow restricting valve in the path of the gas flow.

18. An instrument according to any one of claims 13 to 17, comprising means for detecting the pressure of the purging gas.

10

19. An instrument according to claim 18, comprising means for inhibiting cooling operation of the instrument if the detected pressure is outside a predetermined range of gas pressures.

15

20. An instrument according to any one of the preceding claims, comprising cryogenic fluid conduit connection apparatus having complementary first and second connectors having cooperable locking means operable to lock the first connector to the second connector; in which:

20

the first connector comprises means for activating the locking means to lock the first connector to the second connector; and

the second connector comprises means for releasing the locking means to release the first connector from the second connector.

25

21. An instrument according to claim 20, in which the first and second connectors comprise at least fluid supply and fluid exhaust conduits, corresponding fluid conduits in the first and second connectors being arranged to mate when the first and second connectors are locked together.

30

22. Apparatus according to claim 20 or claim 21, in which the first and second connectors comprise complementary bayonet locking formations.

23. Apparatus according to claim 22, in which the second connector comprises a

WO 95/30379

PCT/GB95/00739

20

sleeve resiliently biased to overlie and prevent release of the bayonet locking formations, the first connector having sleeve displacing formations to displace the sleeve when the first and second connectors are locked together.

5 24. Apparatus according to claim 23, in which the first connector comprises guard formations for at least partially shielding an operator's hand when the first and second connectors are locked together.


10 25. Apparatus according to any one of claims 20 to 24, comprising a heater for heating at least one of the first connector and the second connector.

25. An instrument according to any one of the preceding claims, in which the cryogenic fluid is liquid nitrogen and the gas is nitrogen gas.

15 26. A cryosurgical instrument comprising:
a probe having a cooling tip operable to be cooled by a flow of a cryogenic fluid within the tip;
means for providing a return path for cryogenic fluid from the probe without the fluid passing into the cooling tip; and
20 means for preventing fluid flow via the return path after a predetermined period of operation of the instrument.

WO 95/30379

PCT/GB95/00739

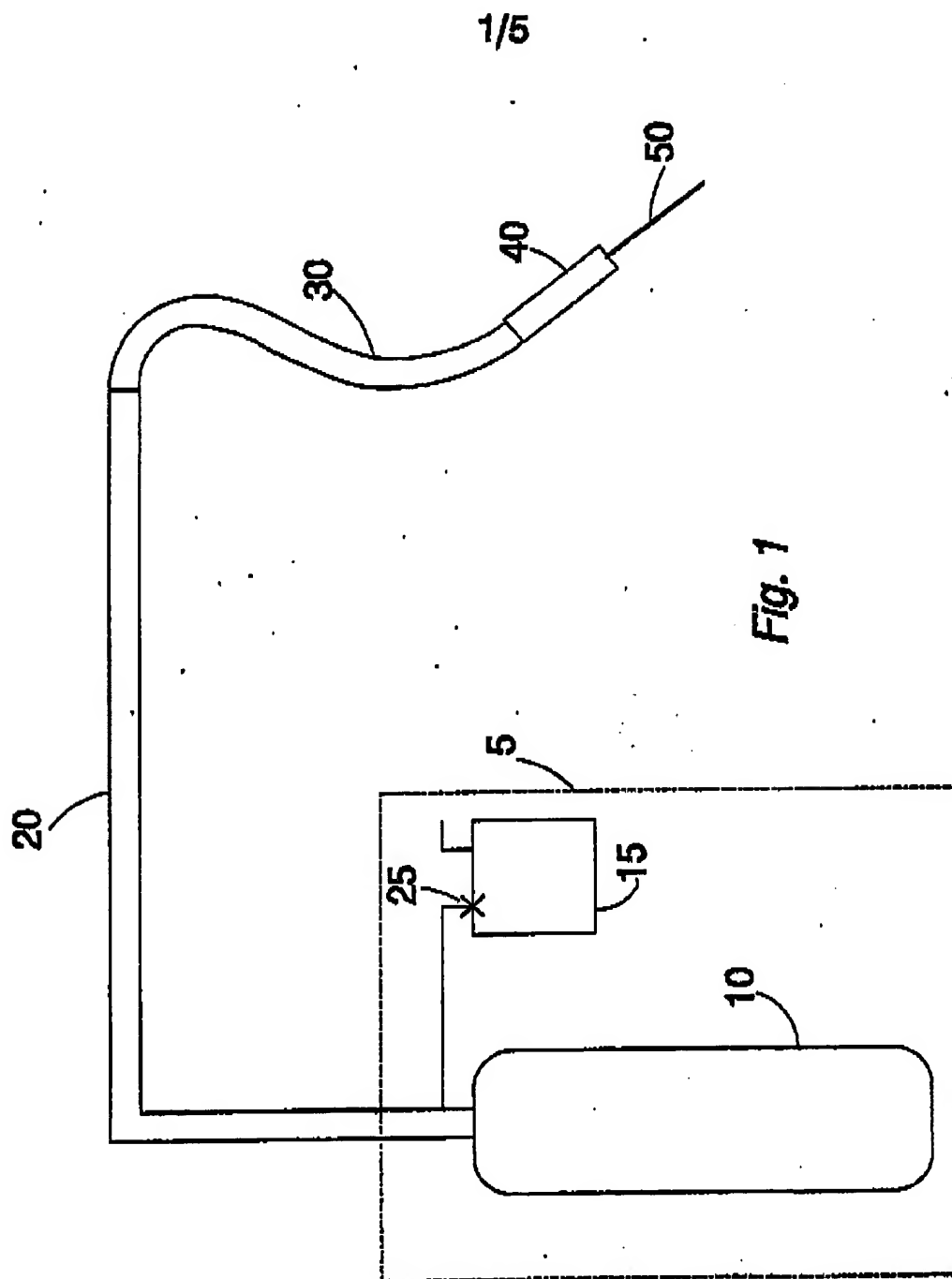

WO 95/30379

PCT/GB95/00739

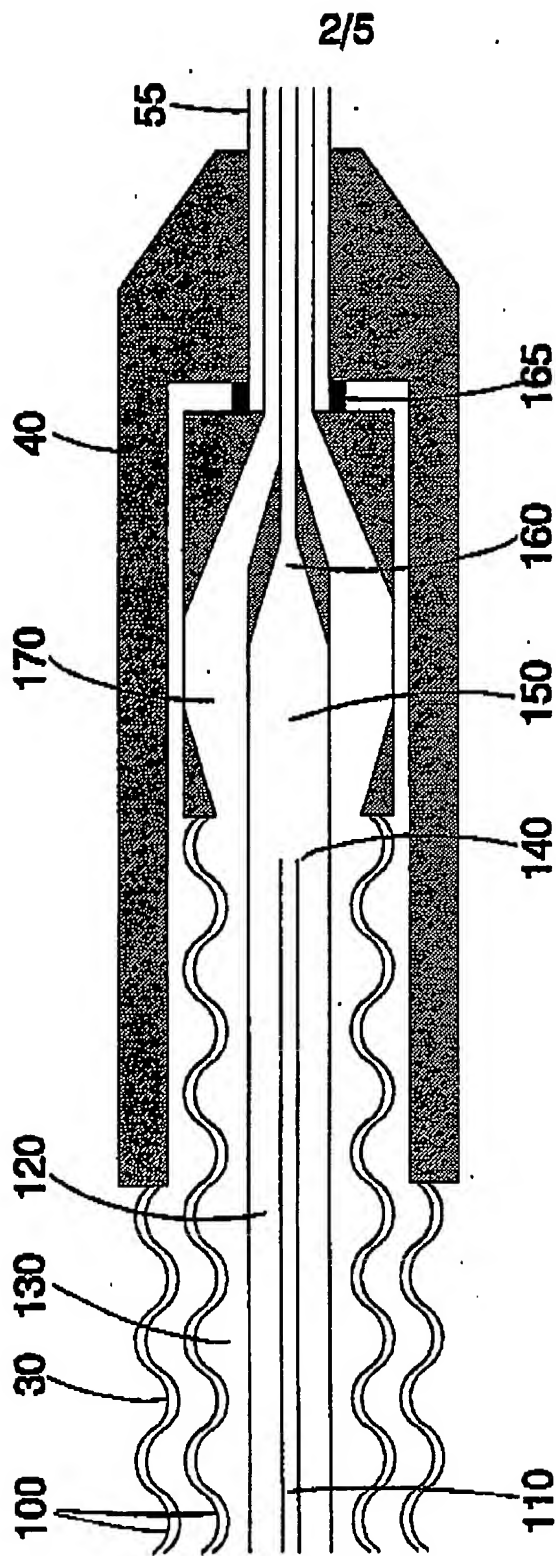


Fig. 2

WO 95/30379

PCT/GB95/00739

3/5

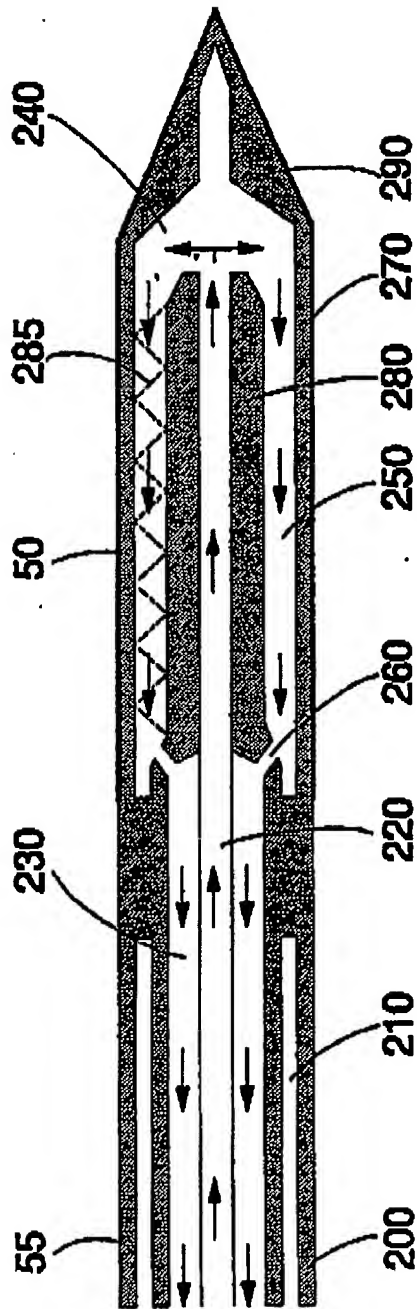


Fig. 3a

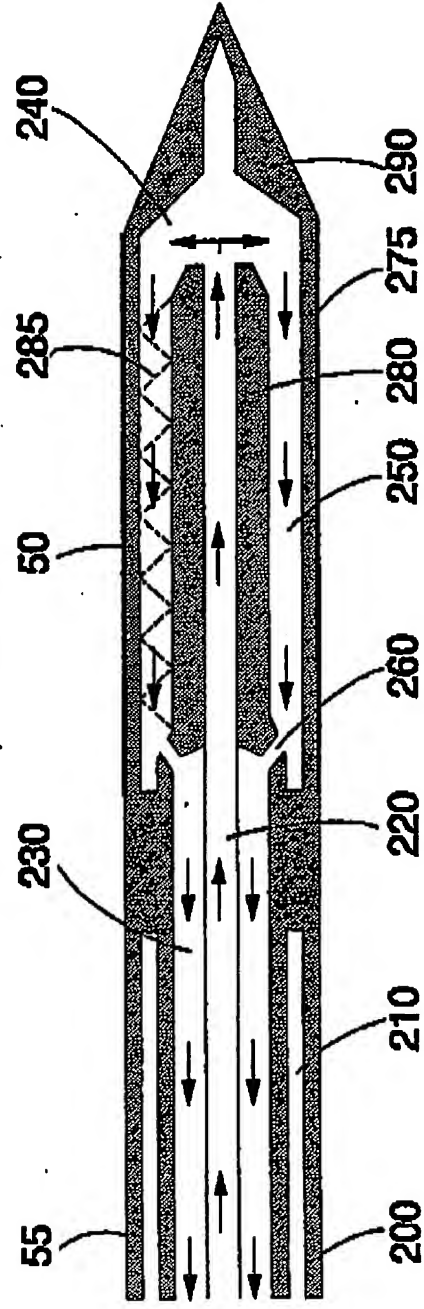


Fig. 3b

WO 95/30379

PCT/GB95/00739

4/5

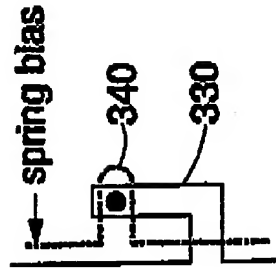


Fig. 5a

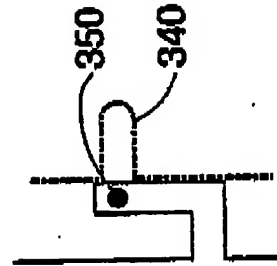


Fig. 5b

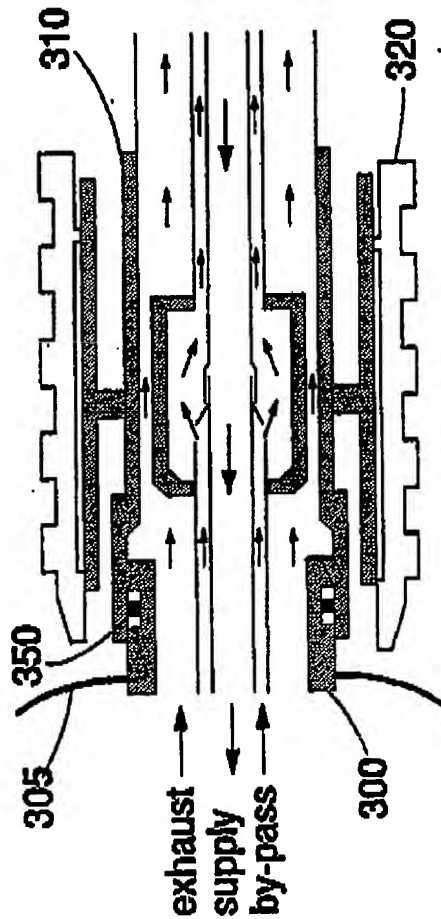


Fig. 4a

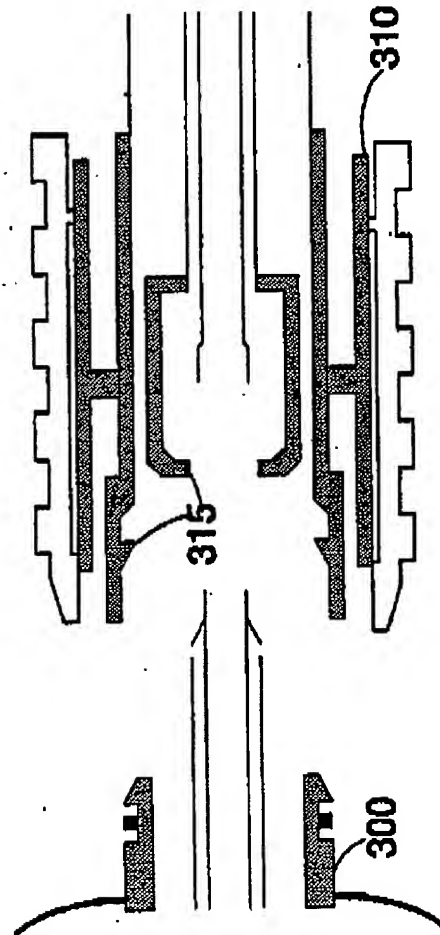



Fig. 4c

Fig. 4b

WFO 95/30379

PCT/GB95/00739

5/5

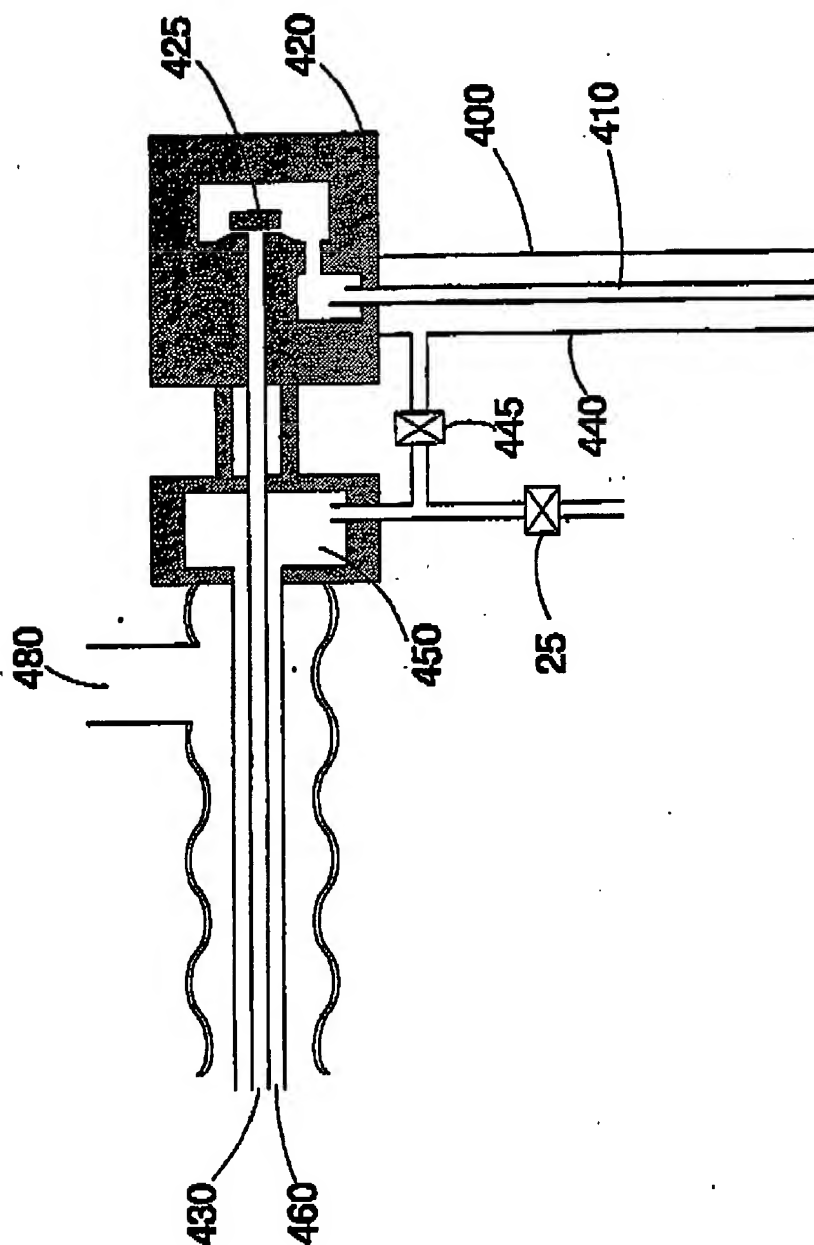

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☐ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.